



# BTL-3000

**BETA**  
**HYDROMASSAGE TUB FOR**  
**LOWER EXTREMITIES**

**USER'S MANUAL**

**CONTENTS:**

<b>1</b>	<b>INSTRUCTIONS FOR USE .....</b>	<b>3</b>
1.1	Outward Appearance .....	3
1.2	Control Elements .....	4
1.3	First Start and Control.....	4
1.3.1	Control of Hydromassage System .....	4
1.3.2	Air Regulation .....	5
1.4	Maintenance and Cleaning .....	5
1.4.1	General Principles of Cleaning.....	5
1.4.2	Cleaning and Disinfection .....	5
1.4.3	Maintenance .....	5
1.5	Principles of Secure Handling.....	6
1.6	Indications .....	7
1.7	Counterindications .....	7
1.8	Terms of Guarantee.....	7
<b>2</b>	<b>TECHNICAL PARAMETERS.....</b>	<b>8</b>
2.1	Transport and Operating Conditions.....	8
2.2	Applicable Standards .....	9

## 1 INSTRUCTIONS FOR USE

The **BETA** series tub is a hydromassage tub for regulable and variable turbulent hydromassage of lower extremities.

The **hydromassage effect** is induced by water flow from large side jets or from microjets which are ergonomically located against the most stressed areas of the upper extremities. Hydromassage is very helpful in treatment of posttraumatic states and injuries. Hydromassage whirlpool bath is an irritation and tonisation procedure and it is also a part of therapeutic physical education. It is effective in elimination of edemas, improvement of mobility of joints and release of contractures. It brings maximum overall relaxation of organism and therefore is optimal for spas, sanatoria, remedial institutions, sporting facilities and resorts of all kinds.

### 1.1 OUTWARD APPEARANCE


The function of the Beta hydromassage tub is based on combination of thermal and mechanical effects of water bath, where the arrangement of microjets ensures turbulent flowing along the whole surface of the lower extremities.

The ergonomic design of the tub enables comfortable position of the patient without useless wetting and with the maximum effect of the bath.

The tub is equipped with a powerful aggregate to which a system of 14 microjets is connected. The tub's system enables pulsation between the inner and outer sides of lower extremities with adjustable interval and smooth air supply of all hydromassage jets. This provides the desired relaxation and/or alleviation of pain in muscles and joints and slows down degenerative processes in joints.

Large enough and clearly organized control panel includes all needed controls which serve for secure and perfect operation of the whole hydromassage equipment. The tub is panelled by covering panels which are adjusted to the shape of the tub. The skeletons of tubs are supplied only in the standard white colour.

## 1.2 CONTROL ELEMENTS

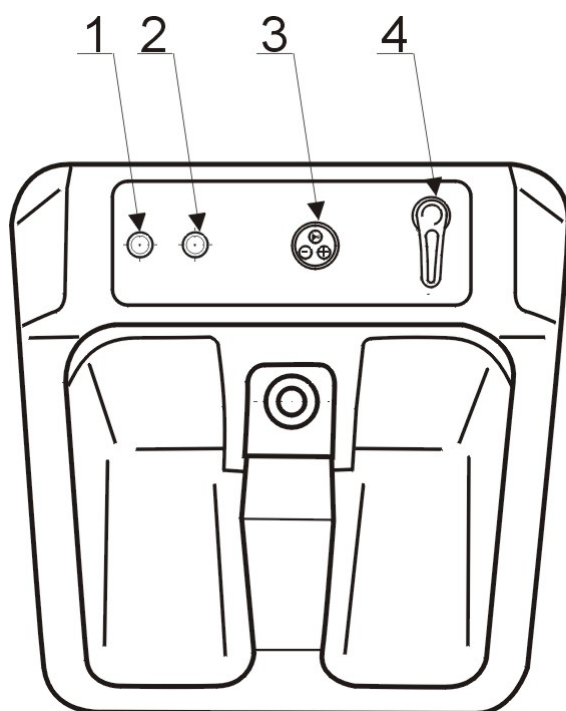

- 1. Air regulation 1
- 2. Air regulation 2
- 3. Electronic control
- 4. Filling of water

Control and operation are intuitive and very comfortable.

- 2x Air regulation – revolving air-supply regulators for outer and inner systems of microjets.
- Hydromassage – electronic control of hydromassage.
- Filling – start of filling of water.

## 1.3 FIRST START AND CONTROL

After successful installation disinfect the tub (see the description below) and rinse its interior.

For filling of water use the lever valve. The tub is filled through the hydromassage system which provides filling in, letting out, overflow and suction of water to the pump.

After filling of the tub (the water top must be at least 2 cm above the topmost microjet) you can start hydromassage.

### 1.3.1 CONTROL OF HYDROMASSAGE SYSTEM

- To start the circulation pump and thus activate hydromassage press the button with the symbol of a jet.
- By the "+" and "-" buttons you can set the interval of switching between the inner and outer sides from 5 to 30 seconds. After each switch-off of the system the time remains saved. If there occurs a power failure or disconnection from the power supply, the system automatically sets the interval to 15 sec.
- Second pressing of the button with the jet symbol switches hydromassage off.

### 1.3.2 AIR REGULATION

The revolving air–supply regulators located on the control panel provide flowing of air into the water systems of hydromassage. Opening of this regulators makes air flow through the jets and microjets together with water. This effect enables to set intensity of water flow to the required level.

### 1.4 MAINTENANCE AND CLEANING

#### 1.4.1 GENERAL PRINCIPLES OF CLEANING

Remove carefully foreign substances (paint, putty, etc.) from the tub's surface, best by a wooden spattle. Never use **steel wool, metallic sponges, knives or scouring washcloth**. Stains of oil or grease can be removed by **denatured alcohol**.

Never use strong solvents such as **acetone, paint thinners, benzine, ammonia or chloride agents, or abrasive polish**.

For metallic parts use special **milk or paste**, like for polishing of a car. Finally wash the cleaned parts by warm water and disinfect them by disinfectants approved by the responsible health officer.

#### 1.4.2 CLEANING AND DISINFECTION

##### Regular cleaning after therapy

- After each therapy let water out, clean the tub by proper cleaner approved by the responsible health officer (e.g. ChiroSan by Bochemie s.r.o.) and rinse it by the shower.
- In the hydromassage system the used water automatically runs out of the hoses (the hoses decline towards the outlet). The outflow from the pump is provided by a drain.
- The equipment's parts which come into touch with the patient should be cleaned by agents approved by the responsible health officer.

##### Regular cleaning and disinfection of pipes 4x per month

Since pollution by mineral sediments and scale considerably affects the activity of internal parts of the system (piping, pump), it is necessary to clean and disinfect the tub regularly.

- a) Fill the tub with warm water, like before bath.
- b) Pour in the water the exact manufacturer-recommended amount of disinfectant (e.g. Savagro by Bochemie s.r.o.).
- c) Run the system for approx. 10 – 15 minutes.
- a) Let water out.
- b) Fill the equipment with warm water and flush the whole system.
- c) Run the system for approx. 10 – 15 minutes.
- d) Let water out.

##### Special cleaning in case of use of non-standard water

If the tub is operated in a place with non-standard water (very hard water, mineral water, etc.) it is necessary to clean it and rinse the hydromassage system at least once a week. The frequency of cleaning depends on quality of water (can be found out by water analysis). For sediments use special agents for dissolution. In case of neglect of proper maintenance and cleaning, complaints will not be accepted.

#### 1.4.3 MAINTENANCE

The service inspection including check of all parameters of the equipment and possible service actions must be performed in intervals complying with the valid law, not longer than 36 months. The inspection is performed by the BTL authorized service department on the basis of the user's order. If the inspection is not done in the stated term the manufacturer does not guarantee the technical parameters and safe operation of the product.

## 1.5 PRINCIPLES OF SECURE HANDLING



This device has applied parts of type B.

The equipment does not use any medicaments or other substances which would be its integral part or would be applied by means of it.



- Before first switch-on of the equipment read carefully the User's Manual.
- The equipment must be professionally installed by an authorized representative of BTL.
- Installation and service instructions are not included in this Manual.
- All staff to use the equipment must be instructed of the way of operation, maintenance and checking of the equipment and of the safety principles.
- The electrical cabling which the equipment will be connected to must be installed and tested according to the existing valid standards (IEC 364). If you are not sure that the mains are completely OK get them inspected by an inspection engineer.
- Check if the parameters of the mains correspond to the requirements of the equipment according to Chapter 2
- The equipment is designed for work in the environment defined in Chapter 2. It must not be used in explosive environment. The equipment must not be used in connection with inflammable anaesthetics or oxidizing fluids (O<sub>2</sub>, N<sub>2</sub>O, etc.).
- Inspect the equipment thoroughly before each use (surface of the tub, functions of displays and controls, etc.); in case of any inconsistency stop using the equipment and contact the authorised service department. If the equipment's behaviour differs from the function described in this Manual stop using the equipment and contact the authorised BTL service department.
- If the equipment shows any defect or if you have doubts about its correct function, terminate the therapy immediately. If you do not determine the source of uncertainty after thorough study of the Manual, contact the authorised service department. If the equipment is used out of accord with this Manual or is used even if it shows functional differences from this Manual, the user is responsible for the damages caused by the equipment!
- Do not dismantle the equipment in any case, removal of protective covers implies the danger of electrical injury.
- All material and parts which come into direct contact with the patient's body must comply with the respective standards related to irritability, allergization, toxicity, genotoxicity and carcinogeneity (ISO 10993-1, ISO 10993-3, ISO 10993-5). The user is responsible for all these materials and parts if not supplied by the BTL equipment supplier.
- The equipment does not use or produce any toxic substances during its operation, storage or transport under the stated conditions.
- Before start of therapy check if all set parameters correspond to your intents.
- To terminate therapy press the respective control element, not the mains switch.
- The equipment and the accessories must not be used in a way out of accord with this User's Manual.
- At work with the equipment use the recommended protective tools.
- The equipment must be placed out of reach of children.
- Main switch use for:
  - switching on and off in operational breaks
  - at repairs and service
  - at weekend downtimes
  - in case of need of fast shut-down
- Do not add to the bath any liquid agents or powders, especially soaps, foams and oils, if not particularly designed for hydrotherapy systems.
- Do not leave persons with restriction in movements, mentally affected persons and children, unattended.
- Near the tub it is forbidden to use any portable electric device! Other electric devices in the room and their parts under voltage must be located and fixed so that they cannot fall into the bath!
- Put the tub into operation before the patient's entry – to prevent him/her from an unpleasant feeling at the first emission of water and air from the jets.
- Do not start the equipment if the tub is not full of water, otherwise the water pump could be damaged.
- After filling up the tub all jets must be under water.
- Before starting of the motor check if there are not any undesired things such as cloth parts, sponges and anything that could be sucked into the pump, which would cause clogging of the pipe and reduction of power or even damage of the pump.
- The equipment contains components which could cause electromagnetic interference.
- It is recommended to separate the patients' rooms from the staff rooms so that the noise level in the staff rooms is reduced (unlike the patients, the staff is exposed to noise for approx. 8 hours a day). In addition it is suitable to

divide the room (at least by curtains) to separate parts, one therapy and one patient each. In case of need extend the anti-noise measures.

- If after many years of operation it is necessary to discard the hydrotherapy equipment it is necessary to contact a specialized company dealing with this activity, or the supplier or manufacturer who will advise you on the process of liquidation, or to discard of it in a way which is usual for this type of devices. The equipment does not contain any toxic materials which could harm the environment in case of normal way of liquidation.

## 1.6 INDICATIONS

- Post-operative states, states after injuries of locomotive organs
- Muscular atrophy
- Post-poliomyelitis states
- Peripheral paresis
- Trophic changes in extremities

## 1.7 COUNTERINDICATIONS

The list of counterindications is the list of cases in which the manufacturer does not recommend application of the selected therapy. Indeed the professional workplaces who are aware of the possible consequences do not need to observe these listed counterindications – all responsibility for use of the therapy, however, lies with them.

- Active TB
- Allergy to therapeutic salts used in the bath
- Skin wounds (including non-bleeding)
- Eczema

## 1.8 TERMS OF GUARANTEE

The manufacturer provides guarantee 24 months from the date of delivery of the hydromassage tub.

The guarantee expires if the equipment has been used out of accord with this Manual or in case of an unqualified intervention in the equipment. The guarantee does not apply to mechanical damage of the skeleton and panelling of the tub, neither to damage of the pump caused by incorrect handling (operation without water).

In case of any defect always contact the authorized BTL service department.

Installation of tubs must be done by qualified personnel with BTL accreditation. In case of "amateurish" installation the supplier does not guarantee for the installation part of the hydromassage equipment and the defects connected with unprofessional installation. Guarantee does not apply for these defects.

## 2 TECHNICAL PARAMETERS

Type	Beta
Material of the skeleton	Acrylate
Maximum volume of the tub [l]	120
Usable volume of the tub [l]	90
Time of filling [min]	1
Time of drainage [min]	1
Weight without water [kg]	55
Dimensions – l x w x h [mm]	920x920x600
<b>Hydromassage system (whirlpool)</b>	x
Number of water jets	0
Number of microjets	14
Air-supply of hydromassage jets	x
Hydromassage pulsation	-
Power supply	230V50Hz
<b>Power supply</b>	
maximum input	800VA
mains voltage	
~ 198 V to 252 V (230 V nominal), alternating	Yes
frequency 50 Hz	Yes
protection class I (according to IEC 536)	Yes
<b>Internal chemical sources</b>	No
<b>Classification</b>	
applied parts of type	B
class according to MDD 93/42/EEC	I Ib

"x" – standard, "o" – optional, "-" – no

### 2.1 TRANSPORT AND OPERATING CONDITIONS

**Identification** BTL-3000 Series System

**Operating conditions**

ambient temperature	+ 10 °C to + 40 °C
relative humidity	30 % to 75 %
atmospheric pressure	700 hPa to 1060 hPa
position	on legs
type of operation	continuous

**Transport and storage conditions**



ambient temperature	- 10 °C to + 55 °C
relative humidity	25 % to 85 %
atmospheric pressure	650 hPa to 1100 hPa
position	on legs
storage time	max. 1 year
other conditions	transport only in the supplied packaging

## 2.2 APPLICABLE STANDARDS

No.	Name	IEC, EN, ISO, MDD
1	Medical electrical equipment Part 1: General requirements for safety	IEC 601-1
2	Amendments to IEC 601-1	A2, A11, A12
3	Medical electrical equipment Part 1: General requirements for safety 2. Collateral Standard: Electromagnetic compatibility. Requirements and tests	IEC 601-1-2
4	Industrial, scientific and medical (ISM) radio-frequency equipment - Radio disturbance characteristics - Limits and methods of measurement	EN 55011
5	Electromagnetic compatibility (EMC) - Part 4: Testing and measurement techniques - Section 2: Electrostatic discharge immunity test - Basic EMC Publication	IEC 61000-4-2
6	Electromagnetic compatibility (EMC) - Part 4: Testing and measurement techniques - Section 3: Radiated, radio frequency, electromagnetic field immunity test	IEC 61000-4-3
7	Electromagnetic compatibility (EMC) - Part 4: Testing and measurement techniques - Section 4: Electrical fast transients/burst immunity test - Basic EMC Publication	IEC 61000-4-4
8	Electromagnetic compatibility (EMC) - Part 4: Testing and measurement techniques - Section 5: Surge immunity test	IEC 61000-4-5
9	Medical devices - Risk Analysis	EN 1441
10	Biological evaluation of medical devices - Part 1: Evaluation and testing	ISO 10 993-1
11	The Medical Devices Directive 93/42/EEC	MDD 93/42/EEC

**Manufacturer**

MEDICAL TECHNOLOGIES CZ a.s.  
Heleny Malirove 11  
169 00 Praha 6  
Czech Republic

**Representative:**

BTL Industries Limited  
Suite 401 Albany House  
324-326 Regents Street  
London  
W1B 3BL  
United Kingdom

tel: +4420 8371 2267

fax: +4420 8371 2270

email: [btl@cobus.net](mailto:btl@cobus.net)

**Sales Department**

BTL zdravotnická technika, a.s.  
Heleny Malirove 11  
169 00 Praha 6  
Czech Republic  
Tel.: (420) 220 517 012, 220 514 753  
Fax: (420) 220 515 520  
E-mail: [sales@btlnet.com](mailto:sales@btlnet.com)  
<http://www.btlnet.com>

*No part of this manual may be reproduced, saved in a research centre or transferred by any means incl. electronic, mechanic, photographic or other records without previous approval from MEDICAL TECHNOLOGIES CZ a.s.*

*BTL is a registered trademark of MEDICAL TECHNOLOGIES CZ a.s.*

*MEDICAL TECHNOLOGIES CZ a.s. operates a policy of continuous development. Therefore, it reserves the right to make changes and improvements to the Product described in this manual without prior notice.*

*The contents of this document are provided "as is". Except as required by applicable law, no warranties of any kind, either expressed or implied, are made in relation to the accuracy, reliability or contents of this document. MEDICAL TECHNOLOGIES CZ a.s. reserves the right to revise this document or withdraw it at any time without prior notice.*